

Our orthopaedic foot supports are available in a version for men and for women, each with or without a microfibre covering.

The contour is suitable for use in our LucRo protective shoes for people with diabetes or rheumatism.

Insole model	Article number	Size	Cover
	971880-000	Women 35-48	-
	971881-000	Women 35-48	Microfibre, gobi light
	971870-000	Men 35-52	-
	971871-000	Men 35-52	Microfibre, taupe



Specialist retailers

Schein Orthopädie Service KG  
 Hildegardstr. 5  
 42897 Remscheid, Germany  
 Tel. +49 2191 910-201  
 Fax +49 2191 910-102  
 remscheid@schein.de  
 www.schein.de



099993-207

04.0937.04.05.17


## Protection for people with diabetes, rheumatism or sensitive feet



The NovaPED soft orthopaedic foot support is specially developed for people with diabetes or rheumatism.

Pressure points can become a serious health risk for people with diabetes. The shape and material combination ensure even pressure distribution and avoidance of pressure peaks, which means the insole is also suitable for prophylactic use by people with diabetes. Due to the good padding characteristics, feet altered by rheumatism can also be successfully cared for.

The foundation is composed of inherently stable structural material with additional stiffening material for shaping. Padding is provided by a resilient, impact-absorbing material combination, which, depending on the version, is covered with moisture-absorbing microfibres. The anatomical shaping of the last includes a shell-like heel guide, a hollowing of the basis V and a retrocapital lateral support.

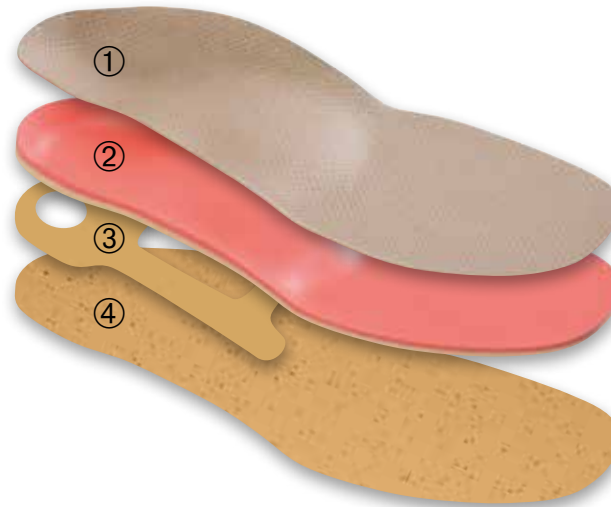

1. **Cover:** Microfibre
2. **Padding layer:** Combination of PE, PU, micro-cork
3. **Stiffening:** THK stiffening
4. **Foundation:** ARU®fant



#### Indications:

- Diabetic foot syndrome (Wagner 0)
- Rheumatic foot

#### Combined with

- Flatfoot
- Splay-foot
- Mildly high arches

